



Antech™ HealthTracks

User Guide

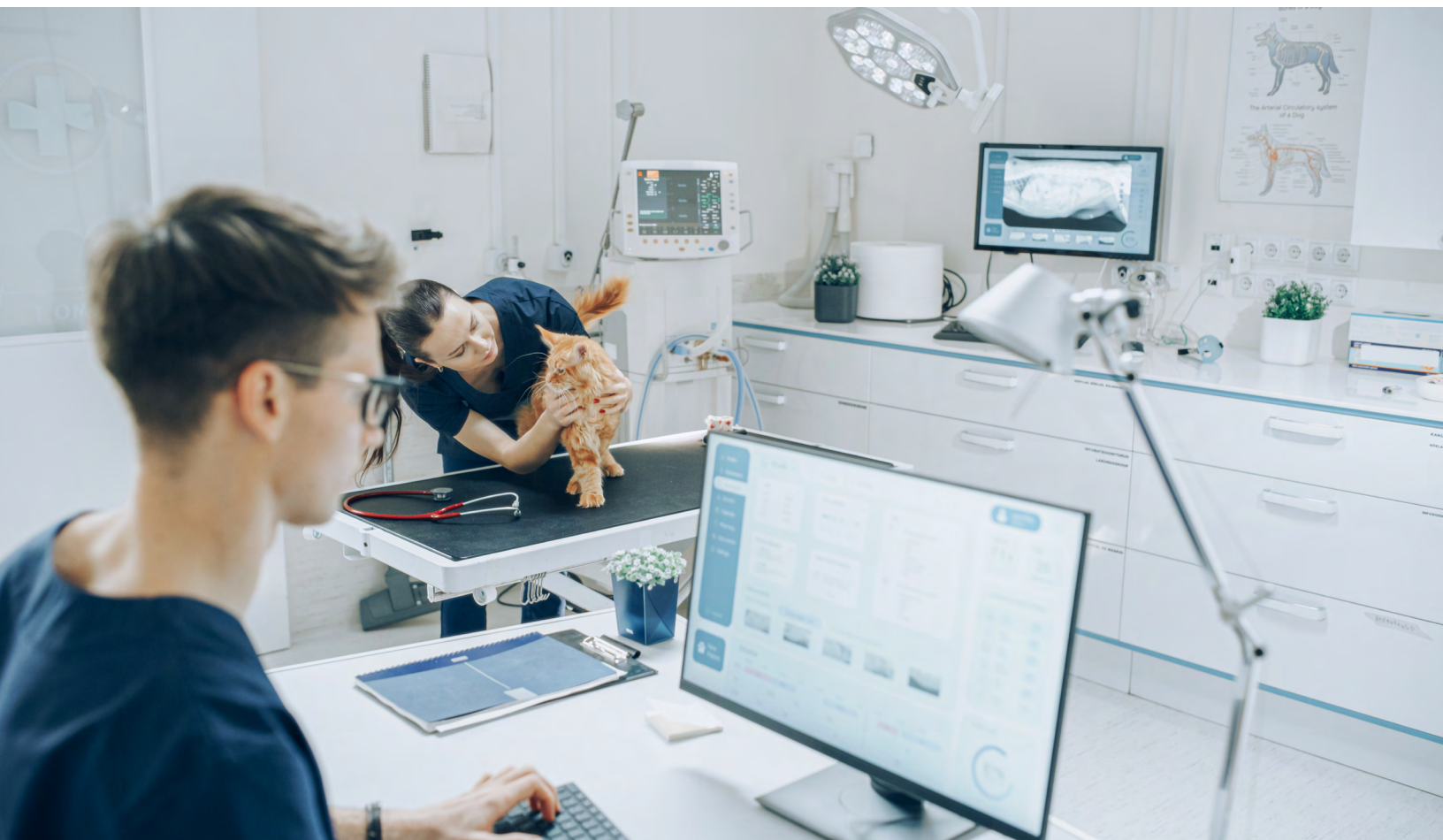


TABLE OF CONTENTS

OVERVIEW 4

ENABLE ANTECH HEALTHTRACKS..... 4

MIGRATION 5

ANTECH HEALTHTRACKS SETUP 7

 Settings7

 Species7

 Breeds..... 8

 Genders 8

 Tests 9

 Code List 9

 Linking Products.....10

SETUP PRODUCT ATTACHMENTS10

TRAVEL SHEET SETUP FOR ANTECH HEALTHTRACKS 12

ANTECH HEALTHTRACKS ULR SETUP 13

ANTECH HEALTHTRACKS WORKFLOW 15

 Product Attachment Workflow 15

 Ordering through HealthTracks 16

 Lab Integrations (Travel Sheet) Workflow 18

ANTECH HEALTHTRACKS LAB RESULTS 20

 Lab Integrations 20

 Antech HealthTracks..... 21

 Lab Screen..... 21

 Medical Records..... 22

ANTECH HEALTHTRACKS MODULE 23

 Antech HealthTracks Module 23

 Customize the Ribbon..... 24

TABLE OF CONTENTS

SECURE ANTECH HEALTHTRACKS 25

Secured Items..... 25


Security..... 26

Overview

Antech has launched a bi-directional integration, now available to all users, enabling immediate order placement and visibility within Antech labs. This integration, known as Antech HealthTracks, includes an embedded portal for streamlined lab order submission and result viewing.

Enable Antech HealthTracks

1. To begin, click  > **Setup** > **Lab Integrations** > **Antech HealthTracks**.
2. Check the option to **Enable Antech HealthTracks**.
 - a. If this is the **first time enabling Antech HealthTracks**, refer to [“Migration” on page 5](#).
 - b. If Antech HealthTracks has been disabled then enabled again, you will need to manually enter your Antech HealthTracks **Credentials: Clinic ID, Username, Password, and Antech HealthTracks ID**. Continue, then, to **Step 3**.
3. Click **Validate** to verify the credentials are correct.
4. When the credentials have been validated, click the **Download** button to pull the latest lab lists from Antech (species, breeds, genders, and tests,). This button is available on all tabs.

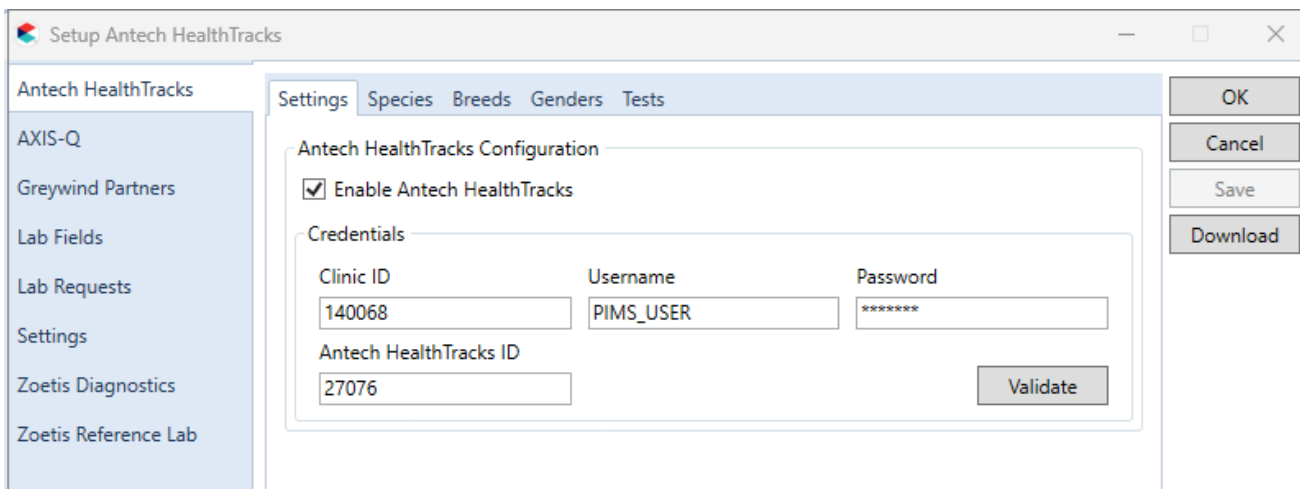


- If the option is not enabled, the Credential fields will be disabled.
- If the credentials are not valid, contact Antech.

Migration


When Antech HealthTracks is enabled, users will be met with a screen to migrate their Antech credentials and settings.

1. To begin, click  > **Setup** > **Lab Integrations** > **Antech HealthTracks**.
2. Check the option to **Enable Antech HealthTracks**. The Migration message appears.

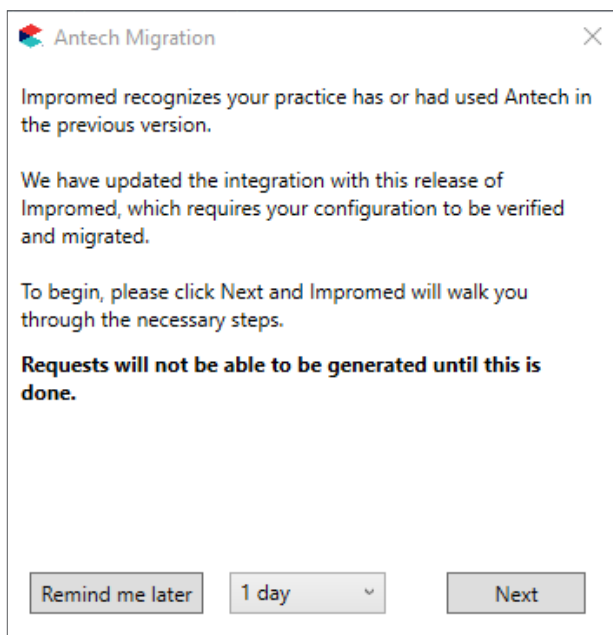

3. Click **Next** to continue with the migration or click Remind me later and set the reminder.
4. After clicking Next, the following screen will appear.

ANTECH HEALTHTRACKS USER GUIDE

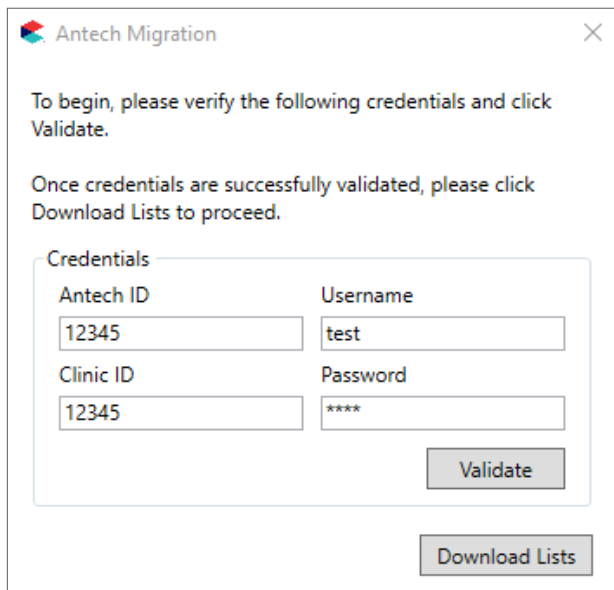

Antech Migration

To begin, please verify the following credentials and click **Validate**.

Once credentials are successfully validated, please click **Download Lists** to proceed.


Credentials

Antech ID	Username
12345	test
Clinic ID	Password
12345	****

Validate

Download Lists

5. Click **Validate** to verify the credentials for your business are correct.
6. When credentials have been validated, click **Download Lists**. A screen will appear with the number of Species, Breeds, and Tests that migrated.

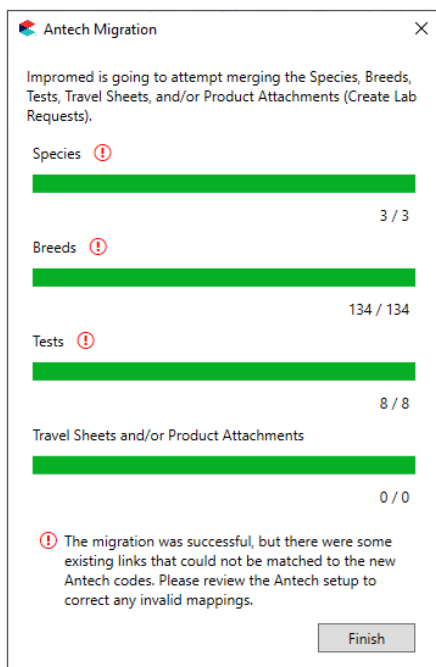

Antech Migration

Impromed is going to attempt merging the Species, Breeds, Tests, Travel Sheets, and/or Product Attachments (Create Lab Requests).

Species ⓘ

3 / 3

Breeds ⓘ

134 / 134

Tests ⓘ

8 / 8

Travel Sheets and/or Product Attachments

0 / 0

ⓘ The migration was successful, but there were some existing links that could not be matched to the new Antech codes. Please review the Antech setup to correct any invalid mappings.

Finish

7. Click **Finish**. The Setup screen will open to verify the mapping is correct, especially for Tests.

Antech HealthTracks Setup

The steps below will describe the setup screens for Antech HealthTracks. The setup page is displayed and saved per business and consists of five tabs: “Settings” | “Species” | “Breeds” | “Genders” | “Tests”. Each tab will be described below.



Impromed will attempt to match items from Species, Breeds, and Genders when the lists are downloaded.


Settings

Impromed will attempt to match items from Species, Breeds, and Genders when the lists are downloaded. However, mapping can be done manually if the lists downloaded mapped incorrectly or did not map at all.

- To begin, click  > **Setup > Lab Integrations > Antech HealthTracks.**

Species


The Species tab is where you map the Impromed species with the Antech HealthTracks species.

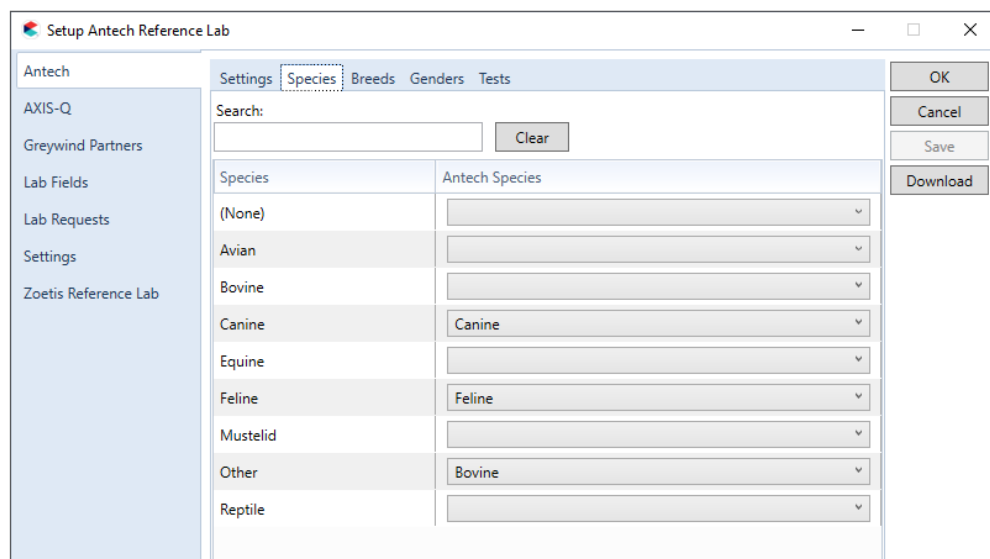

Species	Antech Species
(None)	
Avian	
Bovine	
Canine	Canine
Equine	
Feline	Feline
Mustelid	
Other	Bovine
Reptile	

Breeds


On the Breeds tab, map the Impromed breeds with the Antech HealthTracks breeds.

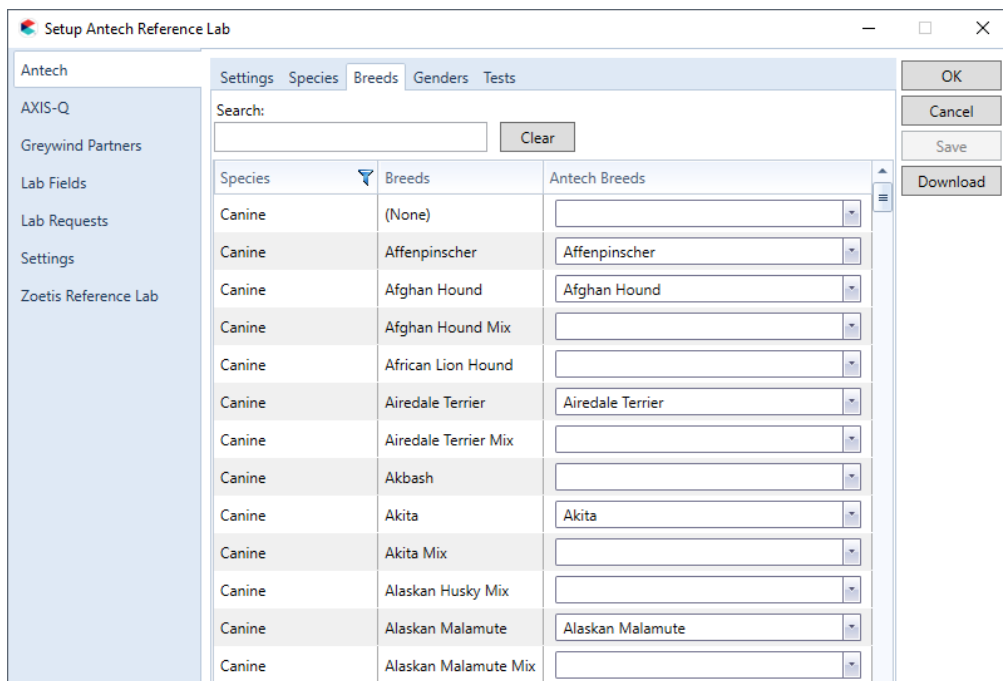

Species	Breeds	Antech Breeds
Canine	(None)	
Canine	Affenpinscher	Affenpinscher
Canine	Afghan Hound	Afghan Hound
Canine	Afghan Hound Mix	
Canine	African Lion Hound	
Canine	Airedale Terrier	Airedale Terrier
Canine	Airedale Terrier Mix	
Canine	Akbash	
Canine	Akita	Akita
Canine	Akita Mix	
Canine	Alaskan Husky Mix	
Canine	Alaskan Malamute	Alaskan Malamute
Canine	Alaskan Malamute Mix	

Genders

The Genders tab is where you map the Impromed genders with the Antech HealthTracks genders.

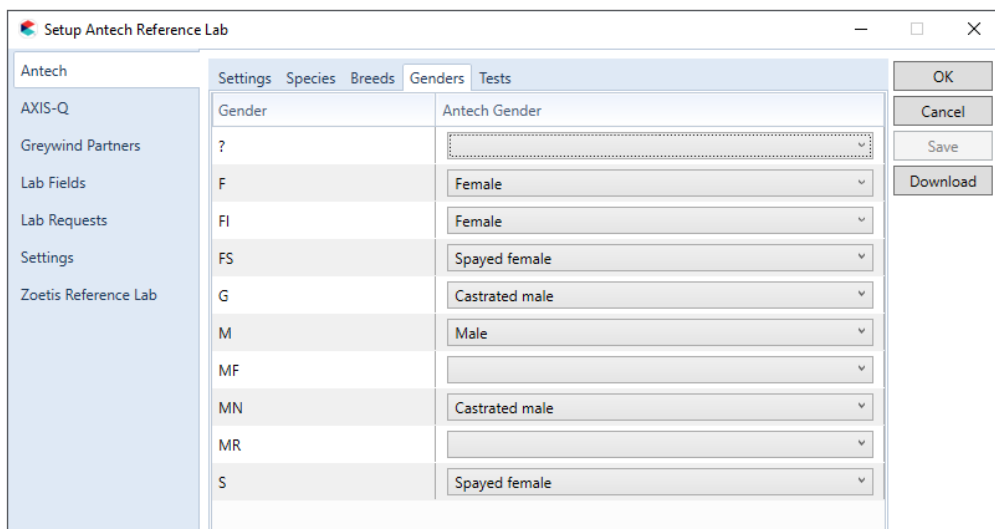


Gender	Antech Gender
?	
F	Female
FI	Female
FS	Spayed female
G	Castrated male
M	Male
MF	
MN	Castrated male
MR	
S	Spayed female

Tests

On the Tests tab, map the Impromed products to the Antech HealthTracks tests. Impromed products mapped to Antech HealthTracks will be the products that appear for selection when creating the Antech HealthTracks travel sheet and product attachments.

Code List

- **Search** - Allows the operator to search by Antech HealthTracks test criteria which filters the table as they type.
- **Lab Test Name** - Test name provided by Antech HealthTracks. You can filter Antech HealthTracks tests to show All tests, Reference only, or In-Clinic.
- **Code** - Code associated to the test provided by Antech HealthTracks.
- **List Price** - The list price associated with the test code provided by Antech HealthTracks.

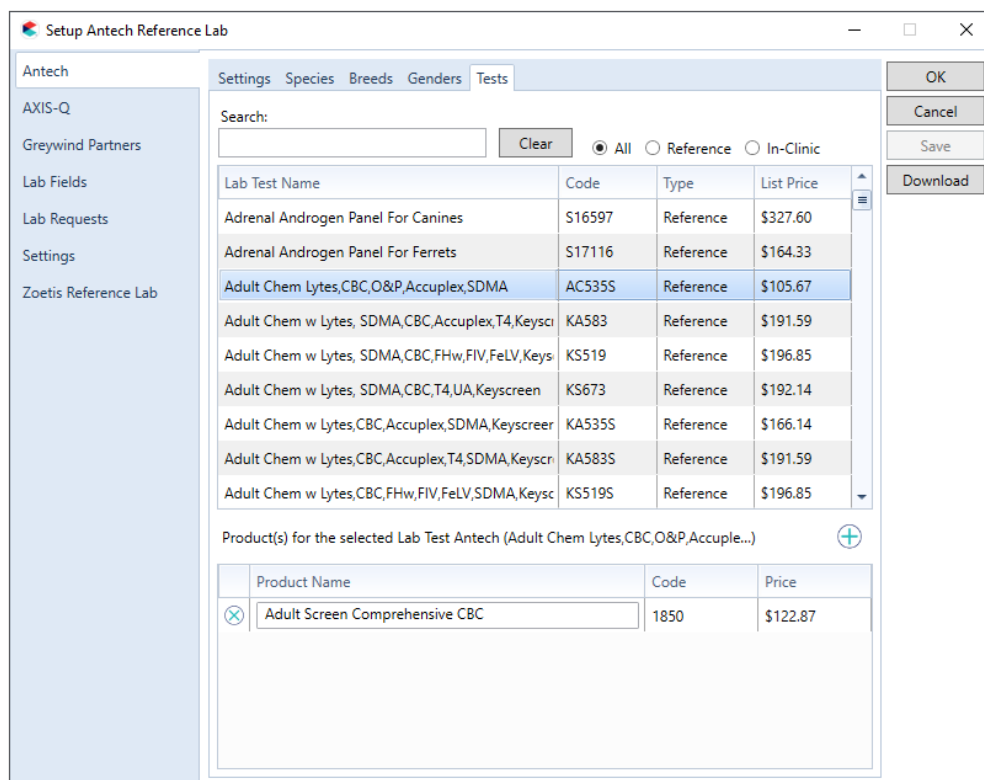

Setup Antech Reference Lab

Antech

AXIS-Q

Greywind Partners

Lab Fields

Lab Requests

Settings

Zoetis Reference Lab

Settings Species Breeds Genders Tests

Search:

Clear

☒ All ☐ Reference ☐ In-Clinic

Lab Test Name	Code	Type	List Price
Adrenal Androgen Panel For Canines	S16597	Reference	\$327.60
Adrenal Androgen Panel For Ferrets	S17116	Reference	\$164.33
Adult Chem Lytes,CBC,O&P,Accuplex,SDMA	AC535S	Reference	\$105.67
Adult Chem w Lytes, SDMA,CBC,Accuplex,T4,Keysci	KA583	Reference	\$191.59
Adult Chem w Lytes, SDMA,CBC,FHw,FIV,FeLV,Keys	KS519	Reference	\$196.85
Adult Chem w Lytes, SDMA,CBC,T4,UA,Keyscreen	KS673	Reference	\$192.14
Adult Chem w Lytes,CBC,Accuplex,SDMA,Keyscreer	KA535S	Reference	\$166.14
Adult Chem w Lytes,CBC,Accuplex,T4,SDMA,Keyscr	KA583S	Reference	\$191.59
Adult Chem w Lytes,CBC,FHw,FIV,FeLV,SDMA,Keysc	KS519S	Reference	\$196.85

Product(s) for the selected Lab Test Antech (Adult Chem Lytes,CBC,O&P,Accuple...)

Product Name	Code	Price
☒ Adult Screen Comprehensive CBC	1850	\$122.87

OK

Cancel

Save

Download

Linking Products


By selecting a test code in the Tests table, the operator is able to use the product grid below the table to add or remove product links.

-  **Add** - Clicking the Add button will add a new row to the bottom of the Product link grid. Search for an Impromed product to link to the selected Antech HealthTracks test.
-  **Remove** - Click to remove a product link/row.
- **Product Name** - This is the name of the product in Impromed.
- **Code** - This is the product code associated to the product in Impromed.
- **Price** - This is the general or base price of the product in Impromed.
 - ◇ Now that the Antech HealthTracks price is displayed, the operator can compare it to the product price and adjust pricing as needed for markup purposes.
 - ◇ Antech HealthTracks pricing would be updated anytime a new list is retrieved.


Setup Product Attachments

Setup products to create an Antech HealthTracks lab request using Product Attachments.

1. Click  > **Setup > Products > Product Attachments**.
2. Choose a product.
3. Click the Attachment Type drop-down menu.
4. Choose **Create Lab Requests**.
5. The request will be made **When the Product is Added to the Invoice**.
6. Click **Maintain**. The Maintain Create Lab Requests for [product] dialog box opens.
7. Select **Antech HealthTracks**.

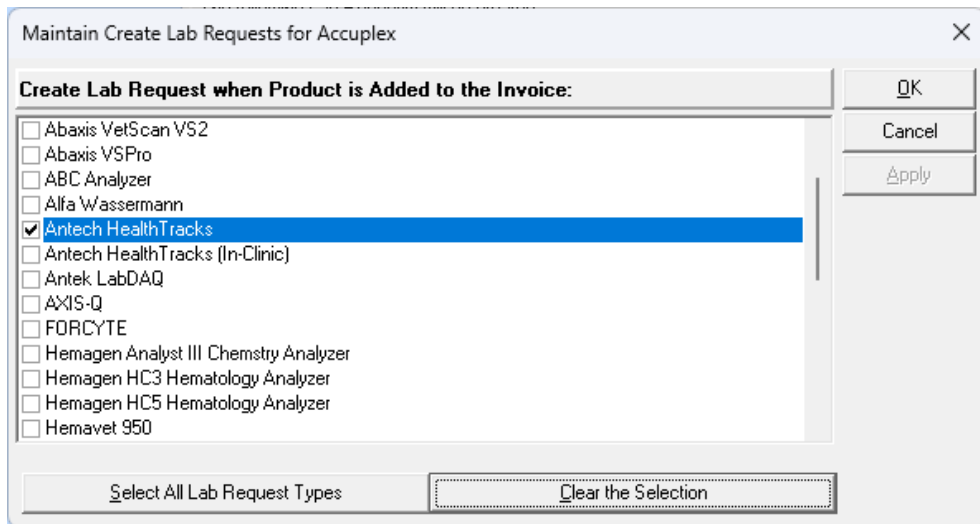

If the product selected is not linked to Antech HealthTracks, the following error will appear:

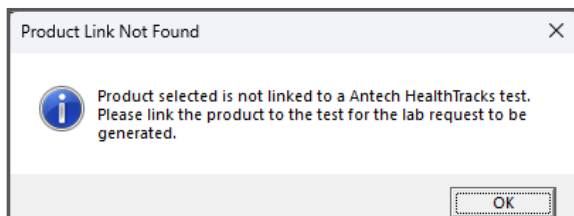


8. Click **OK**.
9. Repeat the steps to attach additional products to the Antech HealthTracks.

Travel Sheet Setup for Antech HealthTracks

For Antech HealthTracks labs, you can create a travel sheet for reference lab requisitions or in-clinic labs. Follow the steps below to create a travel sheet for Antech HealthTracks tests.

1. Click  > **Setup > Travel Sheets**. The Setup Travel Sheets window opens.
2. Click **Create** under **Mode**.
3. Enter the **Name** for the travel sheet.
4. Under **Lab**, select **Antech HealthTracks**.

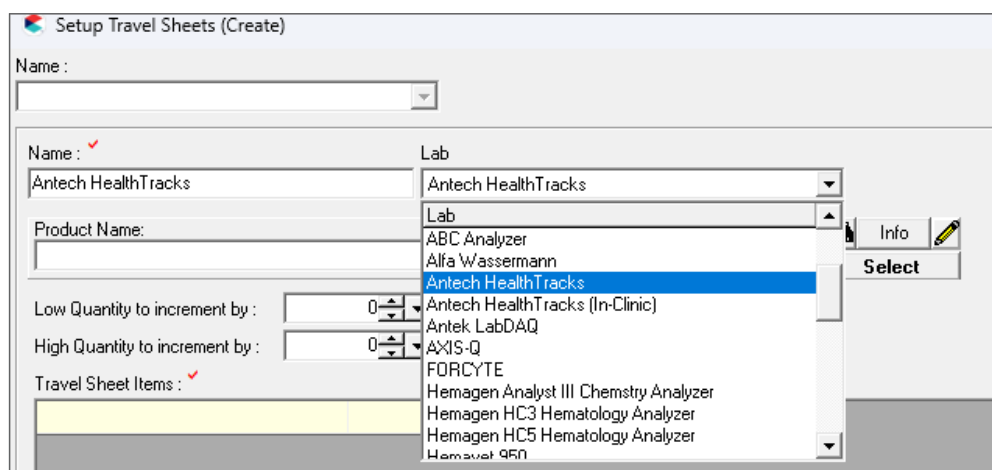

Antech (In-Clinic) requires analyzers to be properly configured by Antech / Heska® for this to work. Until then, In-Clinic is not recommended to be created.

5. From the Product Name, Code, or BarCode, click **Select** for the product to be placed on the travel sheet. **Note:** *The products in the drop-down list are products linked to a lab code from “[Linking Products](#)”.*

Setup Travel Sheets (Create)

Name :

Name : Lab:

Product Name: Code: BarCode:

Name	Code	BarCode	Unit of Measure	Price (qty=1)
Accuplex	1851	1851		\$ 85.00
Adult Screen Comprehensive CBC	1850	1850	Each	\$ 122.87
Biopsy - Muscle or Bone	SM1305	SM1305	Each	\$ 66.13
BUN (Chemistry Analyzer)	L624	L624	Each	\$ 26.45
Canine Diarrhea Panel-Comprehensive	2625	2625	Each	\$ 168.14
CBC & Differential	L628	L628	Each	\$ 51.58
Creatinine	L640	L640	Each	\$ 26.45
E ALT (SGPT)	ET030	ET030	Each	\$ 71.42
E Urine Culture (MIC) Add-On	EADD210	EADD210	Each	\$ 71.42

- Click **Add to Sheet** then repeat the steps to add additional products.
- Click **Apply** to save changes. Click **OK** to save changes and close Setup Travel Sheets.

Antech HealthTracks ULR Setup

Antech HealthTracks is a lab service in the Universal Lab Reader (ULR).

- From the Windows® **Start** menu, go to the **Universal Lab Reader** folder > **Configure ULR**.

Setup Universal Lab Reader v23.11.2.59

File Options

Working Folder (not a Network Path):

SQL Server:

Database Name:

CVET_Improved


Lab Services:

Antech

Running

ANTECH HEALTHTRACKS USER GUIDE

2. Verify the **Working Folder**, **SQL Server**, and **Database Name** have been configured. If not, those fields need to be filled in before proceeding.
3. Click **Add Lab Service**.

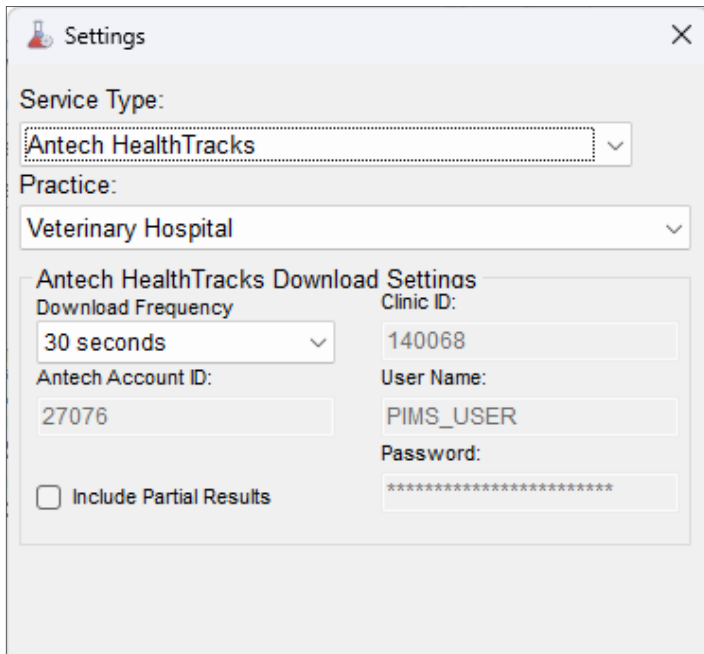

Settings

Service Type:
Antech HealthTracks

Practice:
Veterinary Hospital

Antech HealthTracks Download Settings

Download Frequency: 30 seconds

Antech Account ID: 27076

Clinic ID: 140068

User Name: PIMS_USER

Password: *****

☐ Include Partial Results

4. Select **Antech HealthTracks** from the Service Type drop-down list.
5. Select the **Practice**.
6. Select the **Download Frequency**.
7. The **Antech HealthTracks credentials** will auto-populate based on the credentials entered in Impromed.
8. Check the option to Include Partial Results, if preferred.
9. Click **OK** to save changes.



The Antech HealthTracks service will only be compatible with an Impromed (database) version that has the Antech HealthTracks integration.

If you attempt to add the service to a database that does not contain the appropriate Antech HealthTracks schema in the database, you will be given a message stating that the service cannot be added to the ULR due to incompatibility with the version you are on.

Antech HealthTracks Workflow

Lab ordering can be done from invoicing through Product Attachments or the Travel Sheet. Lab tests can also be ordered directly from the Antech HealthTracks portal within Impromed. All three options are covered below.



A pathology test can be the only test on a lab order for Antech HealthTracks. Operators will need to create a separate order for other tests.

Product Attachment Workflow


If a product has a create lab request attachment for Antech HealthTracks, you are able to order through invoicing. If you have not created product attachments to create lab requests, follow the steps under **“Antech HealthTracks ULR Setup” on page 13.**

1. Click  > **Modules > Invoices.**
2. Select the patient.
3. Add the products with product attachments to create Antech Lab Requests.
4. Conclude the invoice.
5. Antech HealthTracks opens.
 - a. Click to submit samples, if required.
 - b. Scroll down the Cart panel and click on **Place Order.**
 - c. The request will appear which you can download or print.


Ordering through HealthTracks

Lab ordering can be done directly through Antech HealthTracks.

1. Click  > **Modules** > **Lab Integrations**.
2. From the **Requests** tab, click **New** and select **Antech HealthTracks**.

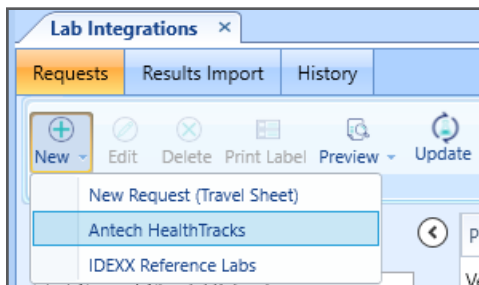

3. The **New Lab Request** appears. Verify the Client, Patient, and Provider are correct.
4. Link the lab to SOAP, Case Number, and/or Problem, if desired.
5. Click **OK**. Antech HealthTracks will open.

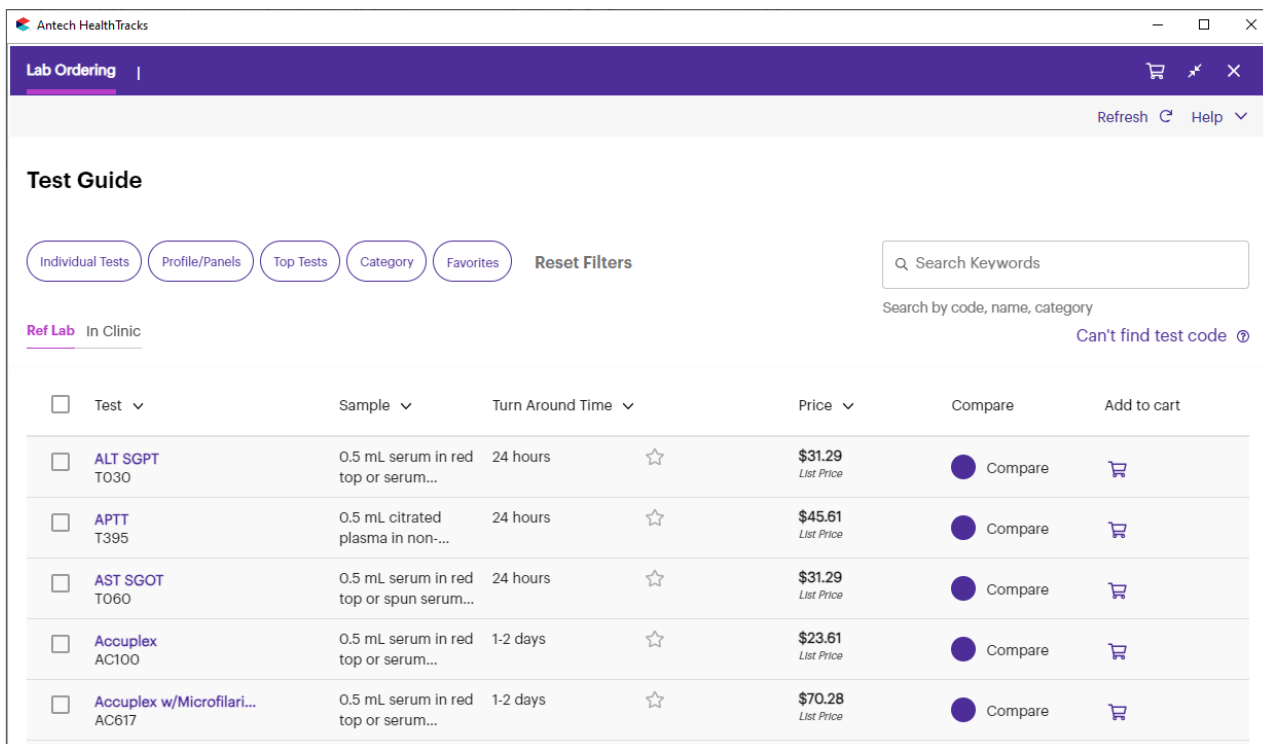

Antech HealthTracks

Lab Ordering | Refresh Help

Test Guide

Individual Tests Profile/Panels Top Tests Category Favorites Reset Filters

Q Search Keywords

Search by code, name, category

Can't find test code ?

☐ Test	Sample	Turn Around Time		Price	Compare	Add to cart
☐ ALT SGPT T030	0.5 mL serum in red top or serum...	24 hours	☆	\$31.29 List Price	☒ Compare	
☐ APTT T395	0.5 mL citrated plasma in non-...	24 hours	☆	\$45.61 List Price	☒ Compare	
☐ AST SGOT T060	0.5 mL serum in red top or spun serum...	24 hours	☆	\$31.29 List Price	☒ Compare	
☐ Accuplex AC100	0.5 mL serum in red top or serum...	1-2 days	☆	\$23.61 List Price	☒ Compare	
☐ Accuplex w/Microfiliari... AC617	0.5 mL serum in red top or serum...	1-2 days	☆	\$70.28 List Price	☒ Compare	

ANTECH HEALTHTRACKS USER GUIDE

6. Locate the test(s) you want to order.

◇ Narrow the search by clicking if the test is a **Ref Lab** or **In Clinic***.

Ref Lab
In Clinic

*In Clinic is only visible if the business uses Antech Analyzers such as Heska.

◇ Use the **Search Keywords** field to help locate tests.

7. Click on the shopping cart icon () to add a test to the Cart.

8. When all tests have been added to the cart, click on the Lab Orders Cart (), located at the top right of the screen, to place the order. The Cart panel opens.

9. If samples are required, click the link to add samples.

Cart

☐ Critical

Account Number

27099

Client ID

--

Doctor

--

Client

--

Pet

--

Pet Age

--

Item Name

Liver Chemistry Panel

List Price

\$67.76

×

Click to add samples*

↗

Total

1 Items

Total Price

\$67.76

Run Priority / Comments

Comments for order (Optional)


10. Scroll down to the bottom of the panel and click **Place Order**.

11. The Antech HealthTracks requisition will appear with the options to close, download, and print.


12. When the order has been placed successfully, the tests are placed on the invoice.

Lab Integrations (Travel Sheet) Workflow

1. Click  > **Modules > Lab Integrations**.
2. From the Lab Integrations module, click the **New** button and choose **New Request (Travel Sheet)**.

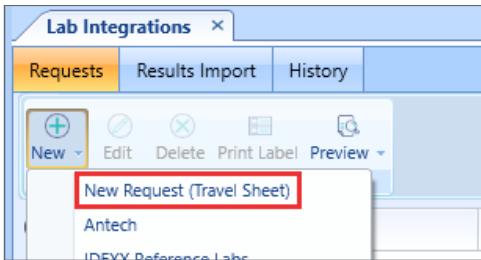


3. On the Lab Requests screen, select the Antech HealthTracks or Antech HealthTracks (In-Clinic) travel sheet you had previously created. If you have not created an Antech HealthTracks travel sheet, refer to [“Travel Sheet Setup for Antech HealthTracks” on page 12](#).

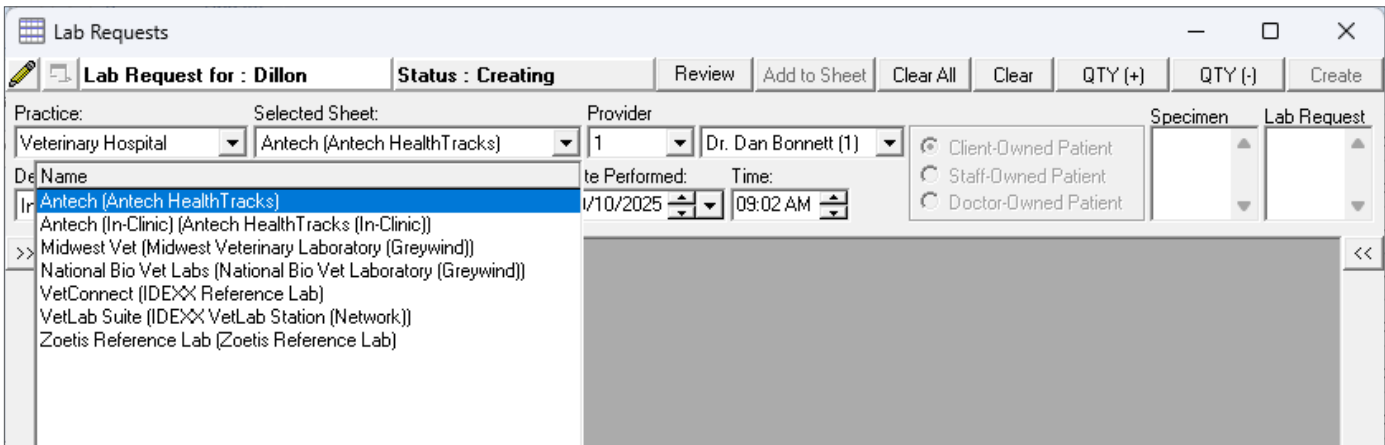

The screenshot shows the 'Lab Requests' window with the status 'Creating'. The 'Lab Request for' field is set to 'Dillon'. The 'Selected Sheet' dropdown is set to 'Antech (Antech HealthTracks)'. The 'Provider' dropdown is set to 'Dr. Dan Bonnett (1)'. The 'Specimen' and 'Lab Request' dropdowns are empty. The 'Patient' dropdown is set to 'Client-Owned Patient'. The 'Date Performed' is set to '1/10/2025' and the 'Time' is set to '09:02 AM'. The 'Name' dropdown is set to 'Antech (Antech HealthTracks)'. The 'Practice' dropdown is set to 'Veterinary Hospital'. The 'Lab Request' list includes: Antech (Antech HealthTracks), Antech (In-Clinic) (Antech HealthTracks (In-Clinic)), Midwest Vet (Midwest Veterinary Laboratory (Greywind)), National Bio Vet Labs (National Bio Vet Laboratory (Greywind)), VetConnect (IDEXX Reference Lab), VetLab Suite (IDEXX VetLab Station (Network)), and Zoetis Reference Lab (Zoetis Reference Lab).

4. Verify the other fields have been filled in correctly.
5. After the order is sent successfully, the order will be presented, allowing the business to download a copy or print the requisition.

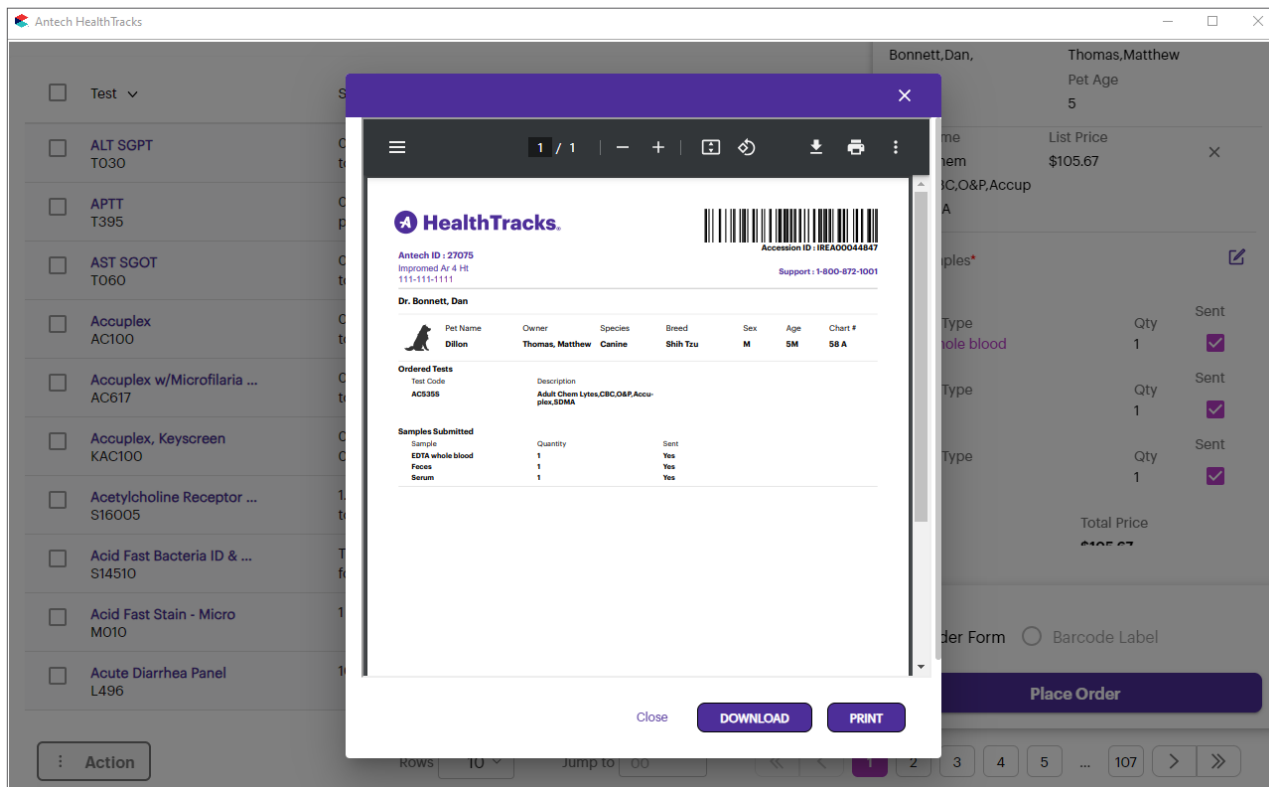

A copy of the order is not received within Impromed. To view the sent order or reprint it, launch the HealthTracks portal > Lab Ordering > click the three dots to the right of the order and choose to download or print the order.

This is where the order can be cancelled in Antech HealthTracks. The order cannot be cancelled from Antech HealthTracks to Impromed. The operator will also need to cancel/delete the order in Impromed.

Antech HealthTracks Lab Results


Imported Antech HealthTracks lab results can be viewed from the Lab Integrations module or from the patient's medical record. Authorized operators can also launch Antech HealthTracks right from Impromed. If the client/patient matches what is in Antech HealthTracks, they will be directed to their latest lab results. If the client/patient does not match, they will be directed to the clinic page.

Lab Integrations

1. Click  > **Modules > Lab Integrations**.
2. To view a lab result, **right-click** on the **imported lab** > **View Result** > choose from **Antech HealthTracks** or **Lab Screen**.




The options to view results via Antech HealthTracks or Lab Screen are also available from the Preview button on the Requests tab.

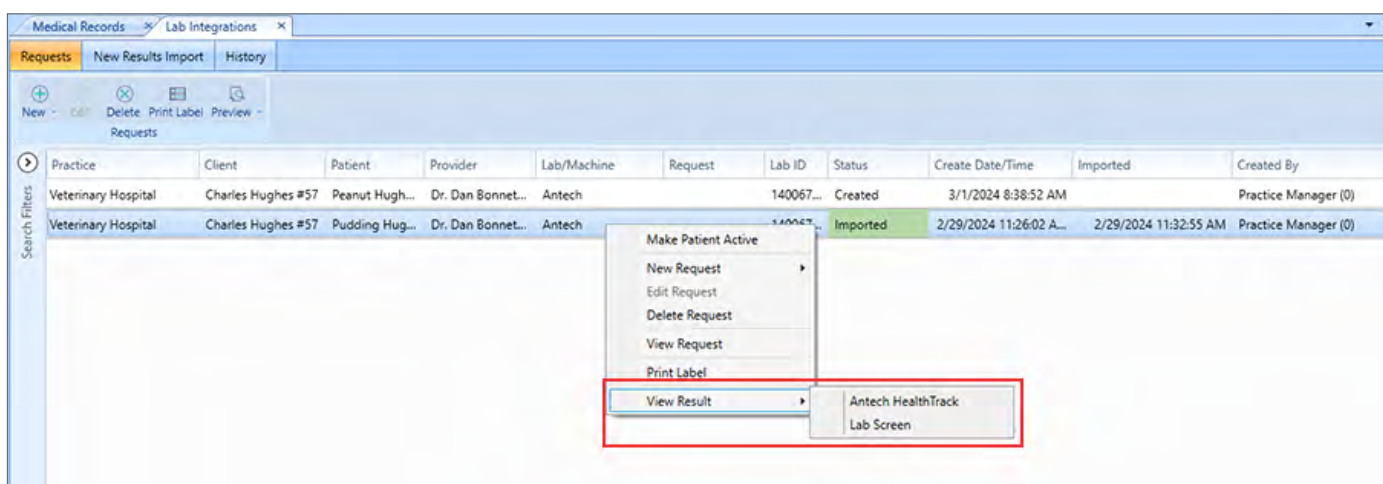

Antech HealthTracks

Antech HealthTracks may provide multiple ways to view the results. In the image below, the view is filtered by Superchem w/SDMA. The option to View All is also available.

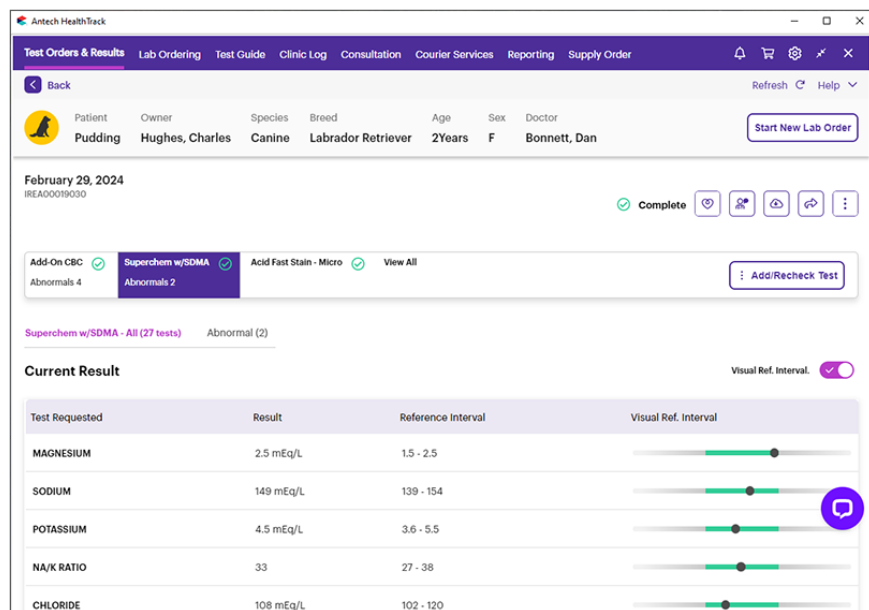


The screenshot shows the Antech HealthTracks web interface. At the top, there's a navigation bar with tabs: Test Orders & Results, Lab Ordering, Test Guide, Clinic Log, Consultation, Courier Services, Reporting, and Supply Order. Below this, a patient summary card displays: Patient (Pudding), Owner (Hughes, Charles), Species (Canine), Breed (Labrador Retriever), Age (2Years), Sex (F), and Doctor (Bonnott, Dan). A 'Start New Lab Order' button is also present. The main section shows a date of February 29, 2024, and a status of 'Complete'. Below this, there are tabs for 'Add-On CBC', 'Superchem w/SDMA', and 'Acid Fast Stain - Micro'. The 'Superchem w/SDMA' tab is selected, showing 'Abnormal (2)' results. A 'View All' button is available. The 'Current Result' section displays a table of test results with visual reference intervals.

Test Requested	Result	Reference Interval	Visual Ref. Interval
MAGNESIUM	2.5 mEq/L	1.5 - 2.5	
SODIUM	149 mEq/L	139 - 154	
POTASSIUM	4.5 mEq/L	3.6 - 5.5	
NA/K RATIO	33	27 - 38	
CHLORIDE	108 mEq/L	102 - 120	

Lab Screen

Below is a sample view of a lab screen.

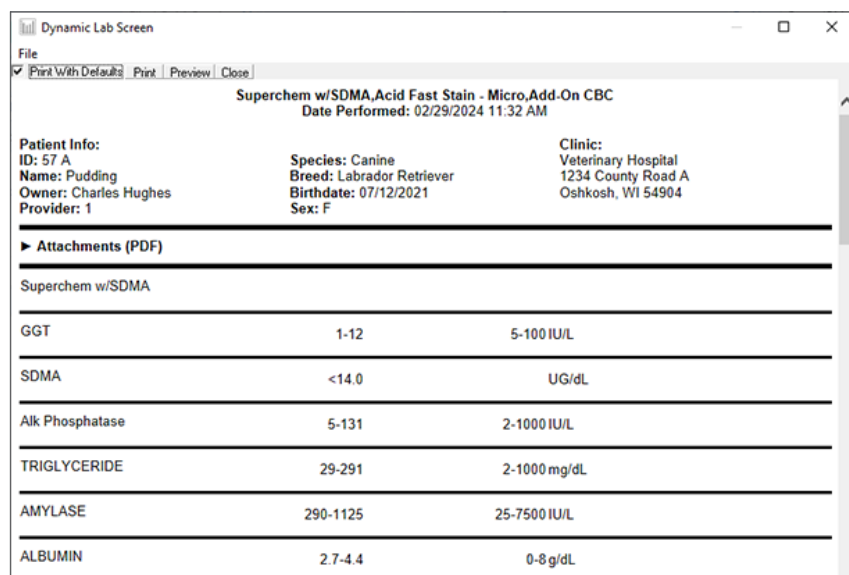

The screenshot shows a 'Dynamic Lab Screen' window. At the top, it says 'Superchem w/SDMA, Acid Fast Stain - Micro, Add-On CBC' and 'Date Performed: 02/29/2024 11:32 AM'. Below this, patient information is displayed: Patient ID: 57 A, Name: Pudding, Owner: Charles Hughes, Provider: 1, Species: Canine, Breed: Labrador Retriever, Birthdate: 07/12/2021, Sex: F, and Clinic: Veterinary Hospital, 1234 County Road A, Oshkosh, WI 54904. There's a section for 'Attachments (PDF)' with a link to 'Superchem w/SDMA'. Below this, a table of lab results is shown.

Test	Result	Reference Interval
GGT	1-12	5-100 IU/L
SDMA	<14.0	UG/dL
Alk Phosphatase	5-131	2-1000 IU/L
TRIGLYCERIDE	29-291	2-1000 mg/dL
AMYLASE	290-1125	25-7500 IU/L
ALBUMIN	2.7-4.4	0-8 g/dL

Medical Records

From a patient's medical record with an imported Antech HealthTracks lab result, simply click on the lab to view the lab screen.

Medical Records										Lab Integrations																			
All	SOAP	Form	Dx	Problem	Image	Lab	Link	Note	Surgery	Invoice	Messages	Dr. Dan	DICOM																
Note	▼		Add Item	Charts	Sheets	Print	DICOM Browser																						
Pudding Hughes: Records 27 of 27										Page 1 of 1										Page 1									
+	02/29/2024	Lab Screen				Superchem w/SDMA,Acid Fast Stain - Micro,Add-On CBC								Provider: 1 Dr. Dan Bonnett DVM															
+		Current Invoice				EM010 E Acid Fast Stain								Provider: 1 Dr. Dan Bonnett DVM															
+		Current Invoice				1850 Adult Screen Comprehensive CBC								Provider: 1 Dr. Dan Bonnett DVM															

To view the results in Antech HealthTracks, click the button on the lab screen.

Medical Records
Lab Integrations

All SOAP Form Dx Problem Image Lab Link Note Surgery Invoice Messages Dr. Dan DICOM

Note Add Item Charts Sheets Print DICOM Browser

Pudding Hughes: Records 27 of 27
Page 1 of 1
Page 1

+ 02/29/2024 Lab Screen Superchem w/SDMA,Acid Fast Stain - Micro,Add-On CBC Provider: 1 Dr. Dan Bonnett DVM

+ Current Invoice EM010 E Acid Fast Stain Provider: 1 Dr. Dan Bonnett DVM

+ Current Invoice 1850 Adult Screen Comprehensive CBC Provider: 1 Dr. Dan Bonnett DVM

Medical Records

All SOAP Form Dx Problem Image Lab Link Note Surgery Invoice Messages Dr. Dan DICOM

Note Add Item Charts Sheets Print DICOM Browser

Vaccine History

Lab Results:

Antech HealthTracks
Preview Email

Superchem w/SDMA,Acid Fast Stain - Micro,Add-On CBC
Date Performed: 02/29/2024 11:32 AM

Patient Info:
ID: 57 A
Name: Pudding
Owner: Charles Hughes
Provider: 1
Species: Canine
Breed: Labrador Retriever
Birthdate: 07/12/2021
Sex: F
Clinic:
Veterinary Hospital
1234 County Road A
Oshkosh, WI 54904

Attachments (PDF)

Superchem w/SDMA

GGT 1-12 5-100 IU/L

SDMA <14.0 UG/dL

Alk Phosphatase 5-131 2-1000 IU/L

TRIGLYCERIDE 29-291 2-1000 mg/dL

AMYLASE 290-1125 25-7500 IU/L


Problems and Diagnosis History Lab Screen

Antech HealthTracks Module

As mentioned previously, Antech HealthTracks can be launched from the Lab Integrations module to create orders and you can launch HealthTracks from Medical Records when viewing lab results. In addition to those areas of Impromed, you can launch Antech HealthTracks as a module and you can add the Antech HealthTracks button to the Ribbon so you can easily launch the portal from within Impromed.

Antech HealthTracks Module

1. Click  > **Modules > Antech HealthTracks**.
2. Antech HealthTracks will open like other Impromed modules.

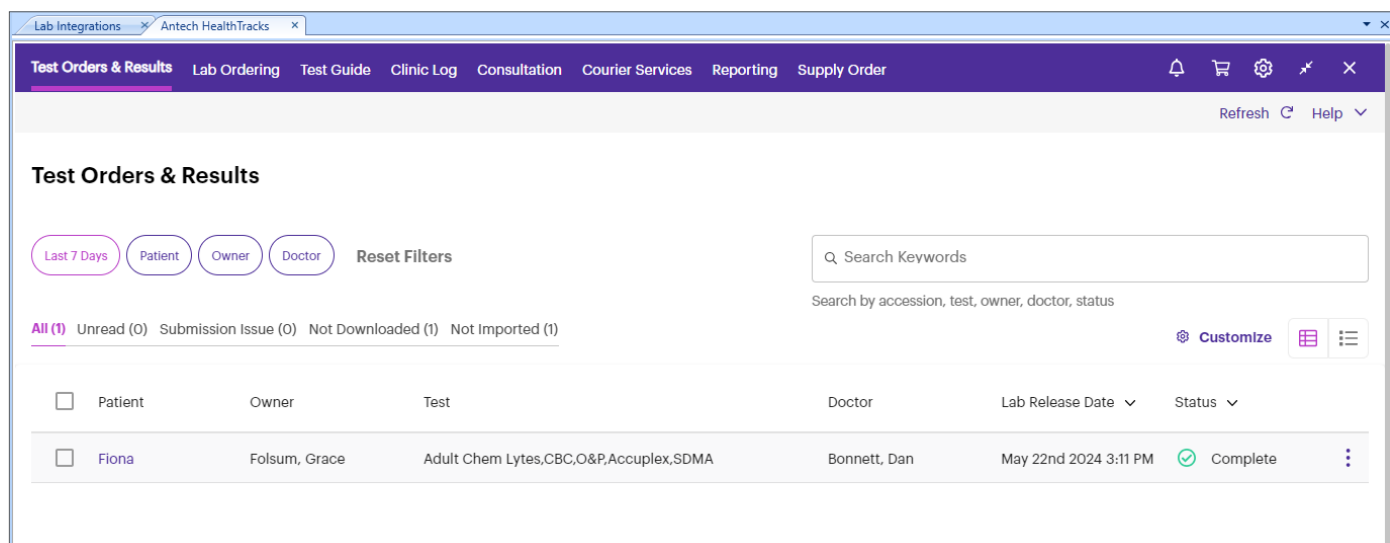

Test Orders & Results

Filters: Last 7 Days Patient Owner Doctor Reset Filters


Search: Search by accession, test, owner, doctor, status

Summary: **All (1)** Unread (0) Submission Issue (0) Not Downloaded (1) Not Imported (1)

☐	Patient	Owner	Test	Doctor	Lab Release Date	Status
☐	Fiona	Folsum, Grace	Adult Chem Lytes,CBC,O&P,Accuplex,SDMA	Bonnett, Dan	May 22nd 2024 3:11 PM	Complete

Customize the Ribbon

1. Click  > **File > Desktop Preferences.**
2. On the Ribbon Bar tab, click on the tab located on the Ribbon to add the Antech HealthTracks button.
3. From the list of Improved Modules, click and drag Antech HealthTracks to where you want the button to appear on the Ribbon. **Note:** In the image shown below, Antech HealthTracks is being added to the Medical Records tab.

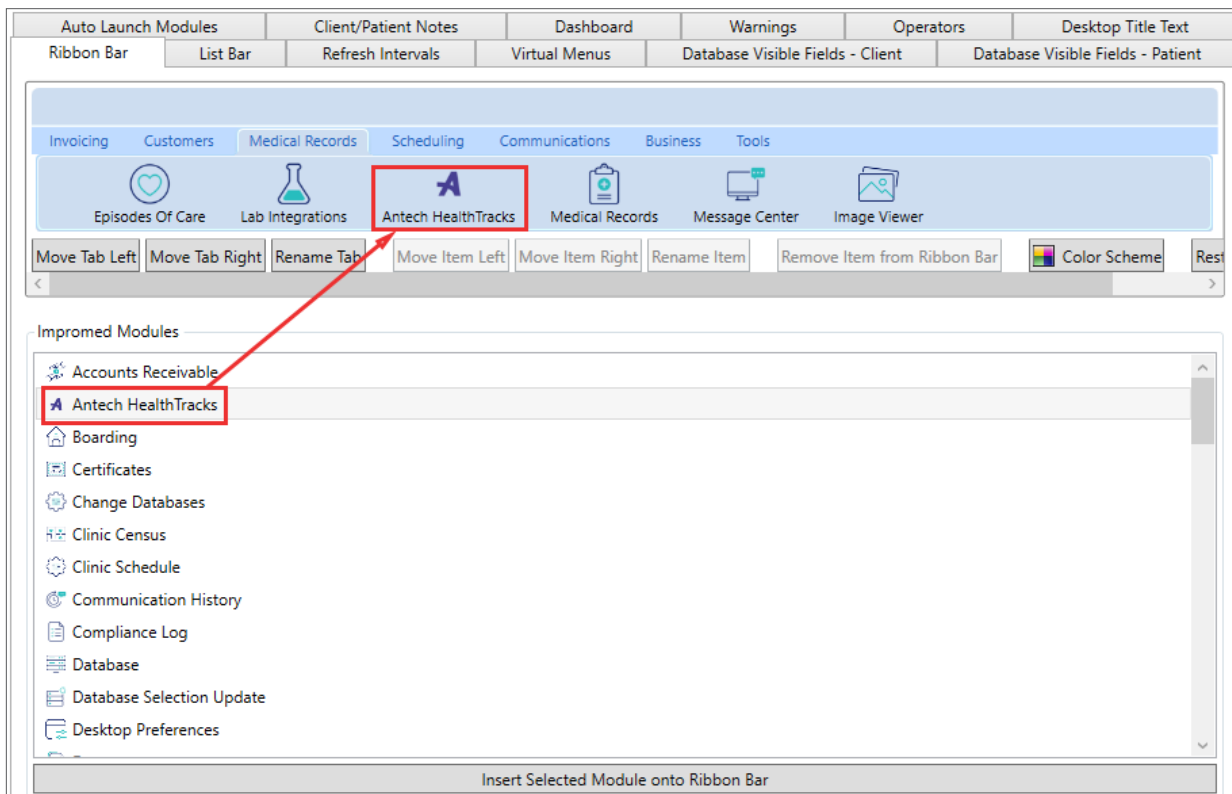


Secure Antech HealthTracks

To restrict employees from opening the Antech HealthTracks, you will need to go into Secured Items and Security.

Secured Items

1. Click  > **Setup** > **Secured Items**.
2. Under Modules, check **Antech HealthTracks** to be able to secure access.

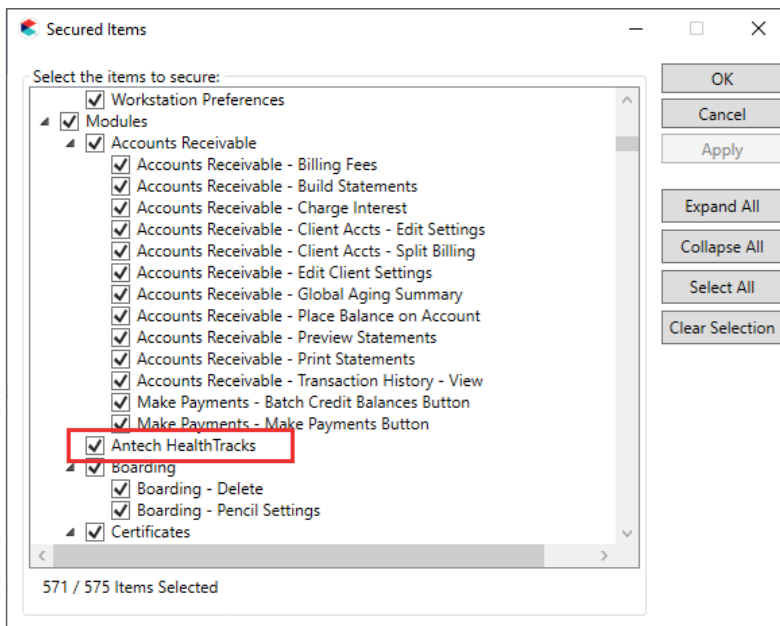


3. Click **OK** to save changes.

Security

Next, you will need to go into Security and make the module secure based on employees selected or employee groups.

4. Click  > **Setup > Security**.
5. Under **Modules**, uncheck **Antech HealthTracks** if the module should not be accessed by the Group Name or Employees selected.

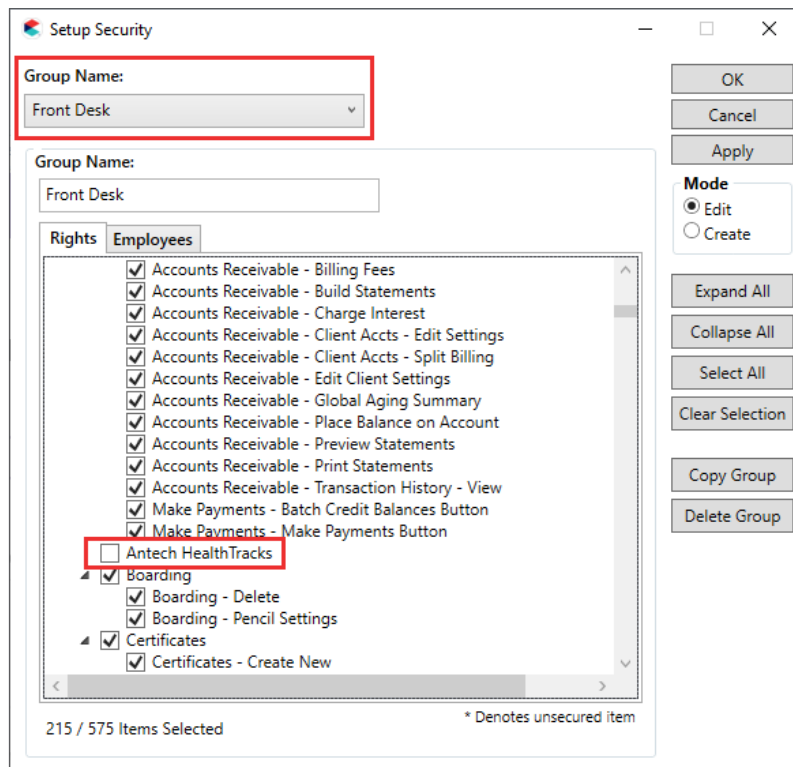

Setup Security

Group Name:
Front Desk

Group Name:
Front Desk

Rights Employees

- ☒ Accounts Receivable - Billing Fees
- ☒ Accounts Receivable - Build Statements
- ☒ Accounts Receivable - Charge Interest
- ☒ Accounts Receivable - Client Accts - Edit Settings
- ☒ Accounts Receivable - Client Accts - Split Billing
- ☒ Accounts Receivable - Edit Client Settings
- ☒ Accounts Receivable - Global Aging Summary
- ☒ Accounts Receivable - Place Balance on Account
- ☒ Accounts Receivable - Preview Statements
- ☒ Accounts Receivable - Print Statements
- ☒ Accounts Receivable - Transaction History - View
- ☒ Make Payments - Batch Credit Balances Button
- ☒ Make Payments - Make Payments Button
- ☐ Antech HealthTracks
- ☒ Boarding
 - ☒ Boarding - Delete
 - ☒ Boarding - Pencil Settings
- ☒ Certificates
 - ☒ Certificates - Create New

215 / 575 Items Selected * Denotes unsecured item

OK Cancel Apply

Mode
☒ Edit
☐ Create

Expand All Collapse All Select All Clear Selection Copy Group Delete Group

6. Click **OK** to save changes.