



Midwest Veterinary Laboratory



IMPROMED®

TABLE OF CONTENTS

- Overview..... 3**

- Setup 3**
 - Settings Tab.....3
 - Species Tab.....4
 - Genders Tab.....5
 - Tests Tab.....6
 - Providers Tab7

- Employee Setup 8**

- Setup Product Attachments 9**

- Travel Sheet Setup 10**

- IULR Configuration 11**
 - Configure the Lab Reader..... 11
 - Add Lab Service..... 12

- Lab Requisition Workflow 13**
 - Product Attachment Workflow..... 13
 - Lab Integrations (Travel Sheet) Workflow..... 13
 - Broken Product Link 14

Overview


Impromed integrates with **Midwest Veterinary Laboratory** to provide the lab workflow from requisitioning to results.

Setup

The steps below will go through the process to enable the use of the Midwest Veterinary Laboratory integration. You will need your credentials for Midwest Veterinary Laboratory.

Settings Tab

1. Click  > **Setup > Lab Integrations > Greywind Partners**.
2. Check the box to **Enable Greywind Partner Integrations**.
3. Check the box to enable the **Midwest Veterinary Laboratory** integration.

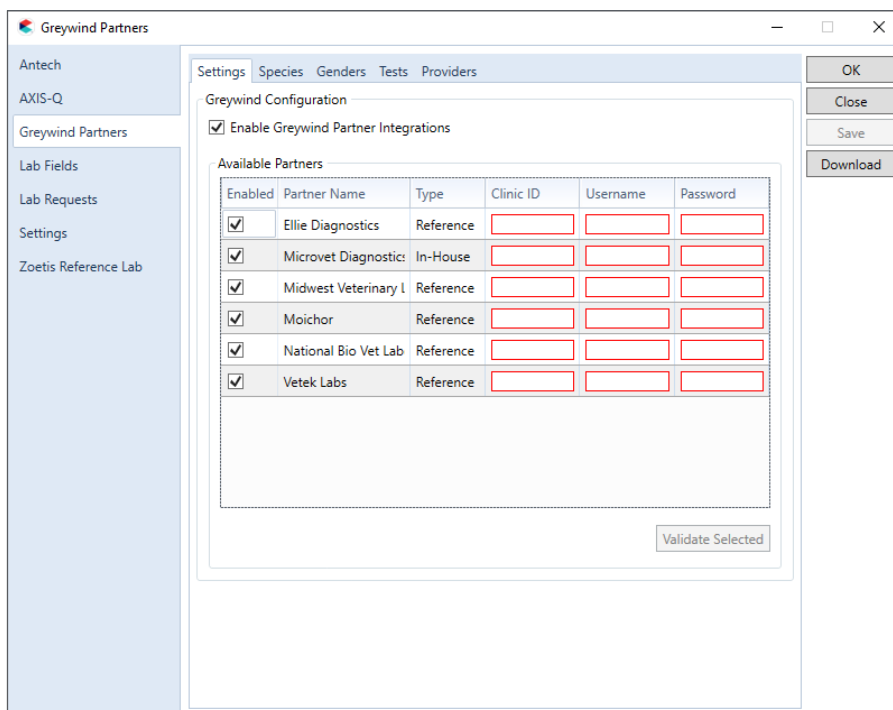

Enabled	Partner Name	Type	Clinic ID	Username	Password
<input checked="" type="checkbox"/>	Ellie Diagnostics	Reference	<input type="text"/>	<input type="text"/>	<input type="text"/>
<input checked="" type="checkbox"/>	Microvet Diagnostic	In-House	<input type="text"/>	<input type="text"/>	<input type="text"/>
<input checked="" type="checkbox"/>	Midwest Veterinary L	Reference	<input type="text"/>	<input type="text"/>	<input type="text"/>
<input checked="" type="checkbox"/>	Moichor	Reference	<input type="text"/>	<input type="text"/>	<input type="text"/>
<input checked="" type="checkbox"/>	National Bio Vet Lab	Reference	<input type="text"/>	<input type="text"/>	<input type="text"/>
<input checked="" type="checkbox"/>	Vetek Labs	Reference	<input type="text"/>	<input type="text"/>	<input type="text"/>

4. Enter the credentials provided to you by Midwest Veterinary Laboratory.
5. Click the **Validate Selected** button. If the credentials entered could not be validated, contact Midwest Veterinary Laboratory to verify your credentials.

6. After validation has been confirmed, click **Download** to map the Species, Genders, Tests available, and Providers.



When Greywind Partner is checked, the Clinic ID, Username, and Password fields will be outlined in red to indicate required fields.

Species Tab

1. On the **Species** tab select **Midwest Veterinary Laboratory** from the Greywind Partners drop-down list to map the Species.
2. After the **Download** button is clicked, the mapping of Species may be done automatically. However, if some species did not get mapped or not mapped correctly, select the species from the Lab Species drop-down list.

The screenshot shows a window titled "Greywind Partners" with a sidebar on the left containing navigation options: Antech, AXIS-Q, Greywind Partners (selected), Lab Fields, Lab Requests, Settings, and Zoetis Reference Lab. The main content area has tabs for Settings, Species (selected), Genders, Tests, and Providers. Under the "Species" tab, there is a "Greywind Partner" dropdown menu set to "Midwest Veterinary Laboratory". Below this is a table with two columns: "Species" and "Lab Species".

Species	Lab Species
Avian	
Bovine	BOVINE
Canine	CANINE
Equine	EQUINE
Feline	FELINE
Mustelid	
Other	
Reptile	

On the right side of the window, there are buttons for "OK", "Close", "Save", and "Download".

Genders Tab

On the Genders tab, map the Improved genders to Midwest Veterinary Laboratory's genders.

1. On the **Genders** tab, select **Midwest Veterinary Laboratory** as the **Greywind Partner**.

If the Download button was clicked, most of the genders should populate. However, this is an opportunity to map genders that weren't downloaded or map those that may have been mapped incorrectly.

2. Click **Save** to save changes.


The screenshot shows the 'Greywind Partners' application window with the 'Genders' tab selected. The interface includes a left sidebar with navigation options: Antech, AXIS-Q, Greywind Partners, Lab Fields, Lab Requests, Settings, and Zoetis Reference Lab. The main area has tabs for Settings, Species, Genders, Tests, and Providers. A dropdown menu for 'Greywind Partner' is open, showing options: Ellie Diagnostics, Microvet Diagnostics, Midwest Veterinary Laboratory, Moichor, National Bio Vet Laboratory, and Vetek Labs. Below this, a table lists various gender codes and their corresponding lab gender mappings:


Code	Lab Gender
FS	
G	
M	Male
MF	
MN	
MR	
S	Spayed

On the right side of the window, there are buttons for OK, Close, Save, and Download.

Tests Tab


To send lab requests to Midwest Veterinary Laboratory, you must map the lab tests from Midwest Veterinary Laboratory to an Impromed product. After doing so, you can add those products to a travel sheet to submit a lab requisition or you can associate the mapped products using Product Attachments to create lab requests when those products are added to an invoice.


1. On the **Tests** tab, select **Midwest Veterinary Laboratory** as the **Greywind Partner**.
2. On the **Lab Test Name** table, select a partner lab test.
3. Below the table, click the **Add** button (). This adds the **Product Name**, **Code**, and **Price** fields to select the Impromed product to match with Midwest Veterinary Laboratory's lab test.

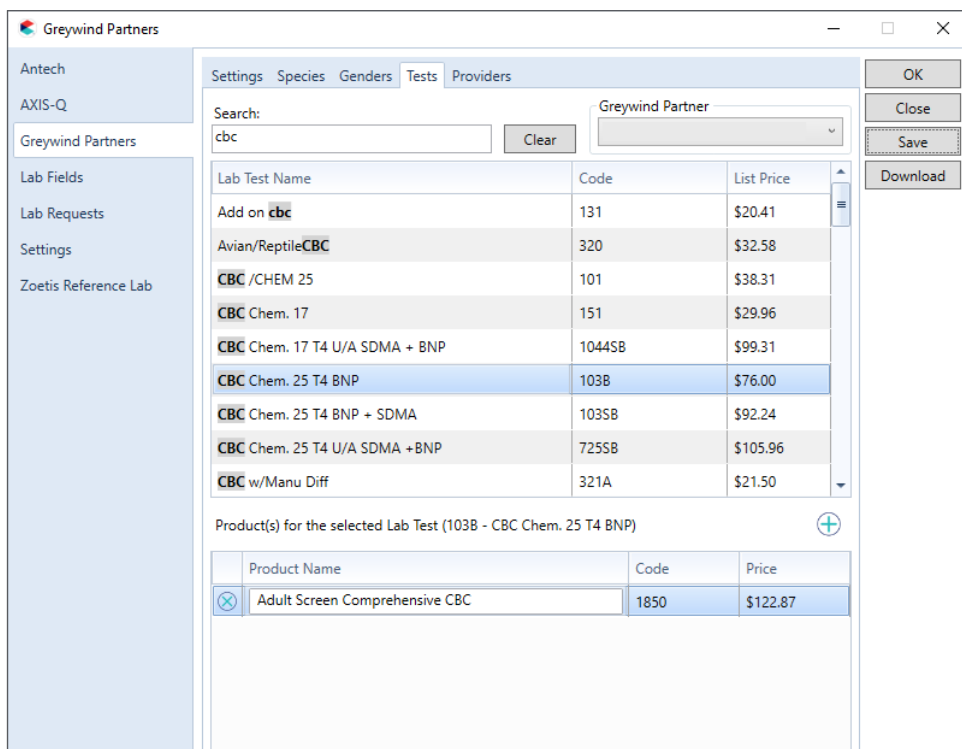

The screenshot shows the 'Greywind Partners' window with the 'Tests' tab selected. The search bar contains 'cbc'. The table below lists various lab tests with their codes and list prices. The 'CBC Chem. 25 T4 BNP' test (code 103B, price \$76.00) is selected. Below the table, a section titled 'Product(s) for the selected Lab Test (103B - CBC Chem. 25 T4 BNP)' contains a table with one row: 'Adult Screen Comprehensive CBC' (code 1850, price \$122.87).

Lab Test Name	Code	List Price
Add on cbc	131	\$20.41
Avian/ReptileCBC	320	\$32.58
CBC /CHEM 25	101	\$38.31
CBC Chem. 17	151	\$29.96
CBC Chem. 17 T4 U/A SDMA + BNP	1044SB	\$99.31
CBC Chem. 25 T4 BNP	103B	\$76.00
CBC Chem. 25 T4 BNP + SDMA	103SB	\$92.24
CBC Chem. 25 T4 U/A SDMA +BNP	725SB	\$105.96
CBC w/Manu Diff	321A	\$21.50

Product Name	Code	Price
 Adult Screen Comprehensive CBC	1850	\$122.87

4. In the **Product Name** field, enter the name of the Impromed product to match Midwest Veterinary Laboratory's test.
5. If an Impromed product was added in error, click the **Remove** button () to remove the product.
6. Click **Save**.

Providers Tab

On the Providers tab, map the Employee / Provider to the Lab Provider for Midwest Veterinary Laboratory. If an employee is not listed as a provider, they may need to be setup as a provider in ["Employee Setup" on page 8](#).

1. On the **Providers** tab, select **Midwest Veterinary Laboratory** as the **Greywind Partner**.

After clicking the **Download** button, the mapping of providers may be done automatically. However, if some providers did not get mapped or not mapped correctly, select the provider from the Lab Provider drop-down list.

2. Under **Lab Provider**, select the employee associated with the **Employee / Provider** in Impromed.
3. Click **Save**.


The screenshot shows the 'Greywind Partners' application window. On the left is a navigation menu with options: Antech, AXIS-Q, Greywind Partners (selected), Lab Fields, Lab Requests, Settings, and Zoetis Reference Lab. The main area has tabs for Settings, Species, Genders, Tests, and Providers (selected). Below the tabs is a search bar with a 'Clear' button and a dropdown menu currently showing 'Greywind Partner'. A table below the search bar has two columns: 'Employee / Provider' and 'Lab Provider'. The table contains four rows of data:


Employee / Provider	Lab Provider
Dr. Dan Bonnett (1)	Dan Bonnett (1400007009)
Dr. J Moore (7)	J Moore (1400007030)
Dr. Lisa Lanzo (2)	Lisa Lanzo (1400007014)
Practice Manager (0)	Practice Manager (1400007029)

On the right side of the interface are buttons for OK, Close, Save, and Download.

Employee Setup

If an employee does not have the rights to submit lab requests, go into **Setup Employees** and give them rights to be a provider of **Lab Screens**.

1. Click  > **Setup** > **Employees** > **Employees**.
2. Select the employee.
3. Click on the **Medical Records** tab.
4. Check **Lab Screens**.
5. Click **Apply**.

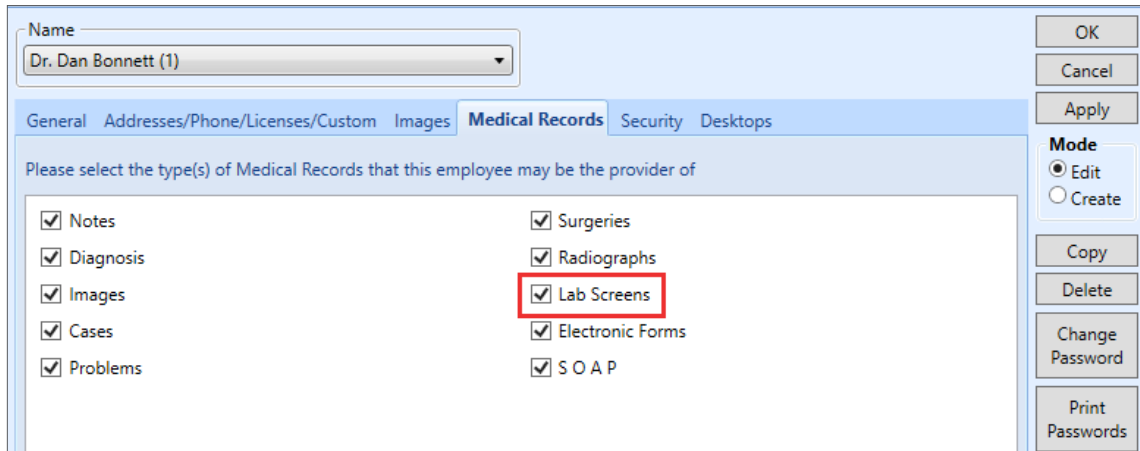


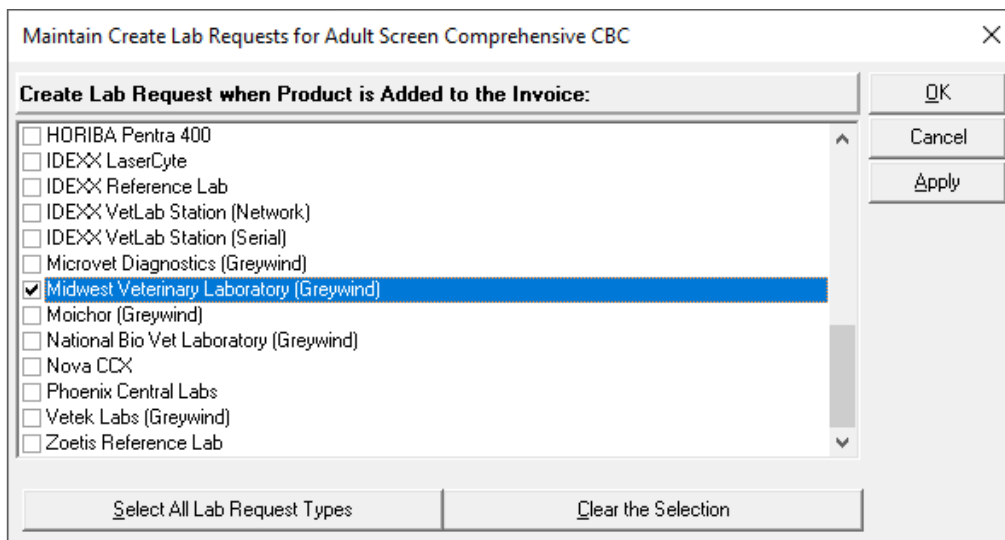
The screenshot shows a software window titled "Employee Setup" for "Dr. Dan Bonnett (1)". The "Medical Records" tab is selected. The window contains a list of medical record types with checkboxes, all of which are checked. The "Lab Screens" checkbox is highlighted with a red rectangular box. On the right side, there are buttons for "OK", "Cancel", "Apply", "Mode" (with "Edit" selected), "Copy", "Delete", "Change Password", and "Print Passwords".

Medical Record Type	Checked
Notes	Yes
Diagnosis	Yes
Images	Yes
Cases	Yes
Problems	Yes
Surgeries	Yes
Radiographs	Yes
Lab Screens	Yes
Electronic Forms	Yes
SOAP	Yes

Setup Product Attachments

Follow the steps below to setup product attachments for Midwest Veterinary Laboratory. Doing so will enable operators to add those products to an invoice which will then generate a lab request to Midwest Veterinary Laboratory.

1. Click  > **Setup** > **Products** > **Product Attachments**.
2. Choose a **Product**.
3. Click the **Attachment Type** drop-down menu.
4. Choose **Create Lab Requests**.
5. The request will be made **When the Product is Added to the Invoice**.
6. Click **Maintain**. The Maintain Create Lab Requests for [product] dialog box opens.
7. Check the box for **Midwest Veterinary Laboratory**.



Maintain Create Lab Requests for Adult Screen Comprehensive CBC

Create Lab Request when Product is Added to the Invoice:

- HORIBA Pentra 400
- IDEXX LaserCyte
- IDEXX Reference Lab
- IDEXX VetLab Station (Network)
- IDEXX VetLab Station (Serial)
- Microvet Diagnostics (Greywind)
- Midwest Veterinary Laboratory (Greywind)
- Moichor (Greywind)
- National Bio Vet Laboratory (Greywind)
- Nova CCX
- Phoenix Central Labs
- Vetek Labs (Greywind)
- Zoetis Reference Lab

Buttons: OK, Cancel, Apply

Buttons: Select All Lab Request Types, Clear the Selection



- If Midwest Veterinary Laboratory is not showing on the Maintain Create Lab list, you need to return to Greywind Setup to enable Midwest Veterinary Laboratory.
- If products were not linked from Midwest Veterinary Laboratory lab tests to Improved products, a warning will display. You will need to return to Greywind Setup - Tests tab to link products.

8. Click **OK**.
9. Repeat the steps to attach additional products to Greywind Partners.


Travel Sheet Setup

This topic presents the steps to create a travel sheet for **Midwest Veterinary Laboratory**.



If you did not link Midwest Veterinary Laboratory lab tests to Impromed products, you will not be able to create a travel sheet.

To create a travel sheet for Midwest Veterinary Laboratory:

1. Click  > **Setup** > **Travel Sheets**. The Setup Travel Sheets window opens.
2. Click **Create** under **Mode**.
3. Enter the **Name** for the travel sheet.
4. Under **Lab**, select **Midwest Veterinary Laboratory** from the drop-down list.



If Midwest Veterinary Laboratory is not showing on the Lab list, you need to return to Greywind Setup to enable the Midwest Veterinary Laboratory.

If products were not mapped from Midwest Veterinary Laboratory lab tests to Impromed products, a warning will display. You will need to return to Greywind Setup > Tests tab to map (link) products.


5. From the **Product Name**, **Code**, or **BarCode**, select the product for the Midwest Veterinary Laboratory lab test.
6. Click **Add to Sheet** then repeat the steps to add additional products.
7. Click **Apply** to save changes. Click **OK** to save changes and close Setup Travel Sheets.

IULR Configuration

To receive lab results from Midwest Veterinary Laboratory, it is necessary to configure the Improved Universal Lab Reader (IULR).


To configure the IULR, follow the steps below. If Midwest Veterinary Laboratory does not appear in the IULR configuration, you will need to enable the partner in Greywind Partners setup.


Configure the Lab Reader

1. Click **Start > [All] Programs > Improved Universal Lab Services > Configure ULR.**
2. Click  for the **Working Folder** location.



The Working Folder must be a local folder and does NOT need to be shared.

3. Select the **C:\Labs** folder.
4. Click **OK**.
5. Enter the **SQL Server** (i.e. pdc1\Impromed - **Note:** this is an Impromed standard naming convention. Your server name may vary).
6. Click  to connect to the SQL Server.
7. Select the **Database Name** (i.e. Impromed).

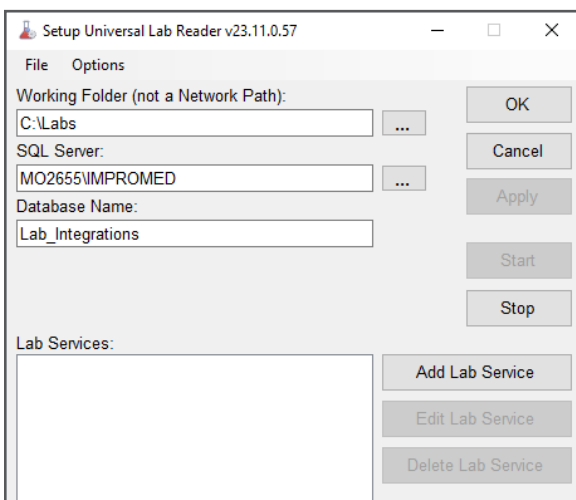


8. Click **Apply**.

Add Lab Service

Once the IULR has been configured, now add the Midwest Veterinary Laboratory lab service.

1. In the Setup Universal Lab Reader, click **Add Lab Service**.
2. From the **Service Type** drop-down list, select **Greywind Partners**.
3. Select the **Practice** from the drop-down list.
4. Configure the Download Settings (**Download Frequency** and **Include Partial Results**).
5. Review the **Enabled Partners**. If a partner is not listed, you need to go to Greywind Partners setup page 3 and enable the partner.

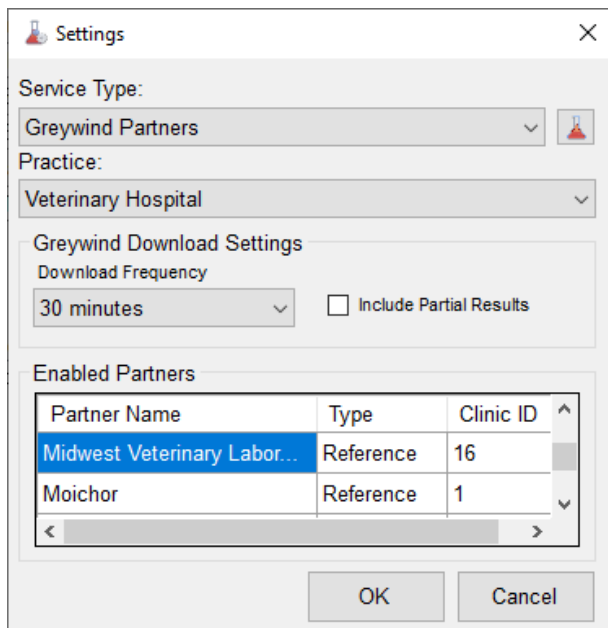

The screenshot shows a 'Settings' dialog box with the following configuration:

- Service Type:** Greywind Partners
- Practice:** Veterinary Hospital
- Greywind Download Settings:**
 - Download Frequency:** 30 minutes
 - Include Partial Results:**
- Enabled Partners:**

Partner Name	Type	Clinic ID
Midwest Veterinary Labor...	Reference	16
Moichor	Reference	1


Buttons: OK, Cancel

Lab Requisition Workflow


This topic covers both ordering labs from the Travel Sheet and Invoicing through Product Attachments.


Product Attachment Workflow

If a product has a create lab request attachment for a Greywind Partner, you are able to order through Invoicing. If you have not created product attachments to create lab requests, follow the steps below the "[Setup Product Attachments](#)" on page 9.

1. Click  > **Modules > Invoices**.
2. Select the patient.
3. Add the products attached to Midwest Veterinary Laboratory Lab Requests.
4. Conclude the invoice.

Lab Integrations (Travel Sheet) Workflow

1. Click  > **Modules > Lab Integrations**.
2. From the Lab Integrations module, click the **New** button, choose **New Request (Travel Sheet)**.

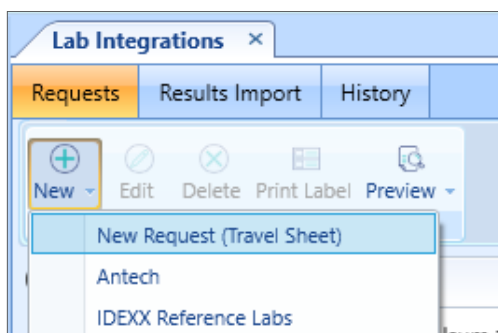

3. On the Lab Requests screen, select the **Midwest Veterinary Laboratory** sheet you created. If you have not created a Midwest Veterinary Laboratory travel sheet, refer to the "[Setup Product Attachments](#)" on page 9.
4. Verify the other fields have been filled in correctly.



If the Provider selected on the travel sheet is not linked to Midwest Veterinary Laboratory, a prompt will appear to select the Provider.

Select Midwest Veterinary Laboratory (Greywind) Provider

Midwest Veterinary Laboratory (Greywind)

Dan Buller (1048)

OK

Cancel

5. After the order is sent successfully, a PDF will be presented, allowing the business to print/save. The request can be viewed by double-clicking on the Requests tab.

Broken Product Link

In the event a product is removed from **Midwest Veterinary Laboratory Tests**, but the product is not removed from the travel sheet, a prompt will display for the operator to select the Midwest Veterinary Laboratory Test to link back to the product when a lab request is being made from the travel sheet.

Select Test for Product

Product to link

Adult Screen Comprehensive CBC

Select test

OK

Cancel

1. Click in the **Select test** field and begin typing the name of Midwest Veterinary Laboratory's test to match the Impromed product.
2. Click **OK** to save changes.

Saving the test here will link the test to the product under Midwest Veterinary Laboratory Greywind Partners Setup.