



# IDEXX® In-House (IVLS)

# User Guide



TABLE OF CONTENTS

Overview..... 3

Advanced Options..... 3

Mapping ..... 4

Select the Device ..... 5

Place the Order ..... 5


Notifications..... 6


## Overview

The new IDEXX IVLS integration has the ability for Avimark to communicate with multiple IVLS devices in the practice. This document will go through the process of configuring the settings.

## Advanced Options

Enter the IDEXX In-House credentials under Advanced Options in Avimark.


1. Click the gear icon (  ) on the speedbar.
2. Search for 'idexx in-house'.
3. Set **Active** to **True**.
4. Enter the **API Password** and **API Username**.
5. Enter the **LIMS ID**.  
The LIMS ID with username and password will let the program know what IVLS machines are in-house.
6. Click **Done**.

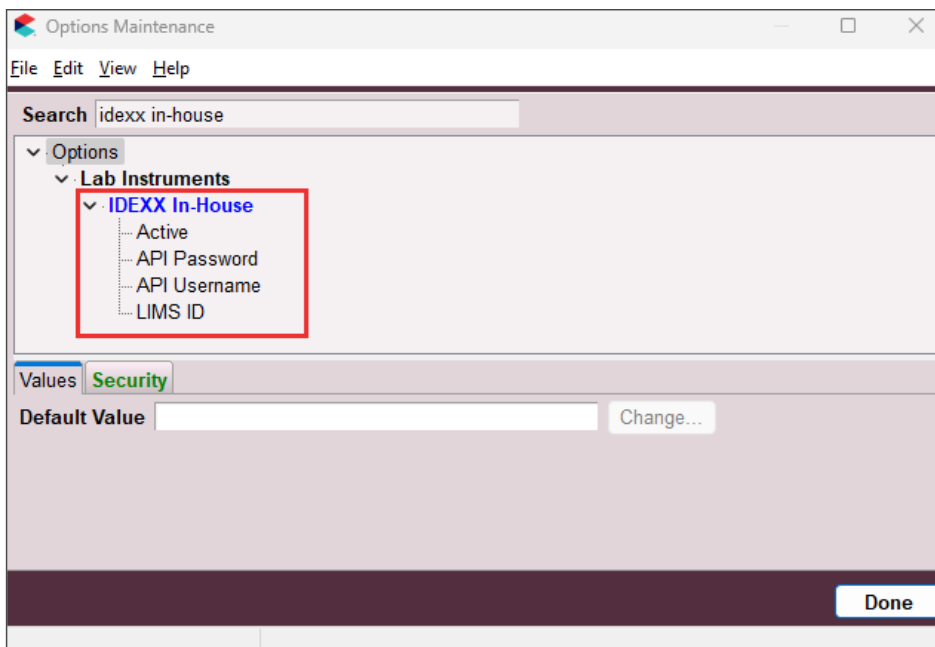

## Mapping

Because the IDEXX IVLS In-House is a new integration, species and breed mapping will be performed automatically. However, mapping can also be done manually, if needed.

1. From the menu, click on **Work with > System Tables**.
2. Select the **Species** or **Breed** table.
3. Double-click on the species or breed that requires mapping.
4. Click on the **Mapping** tab.
5. Under **Code** next to IDEXX In-House, double-click to get the list and select the mapping code for the species or breed selected.
6. Click **OK**.

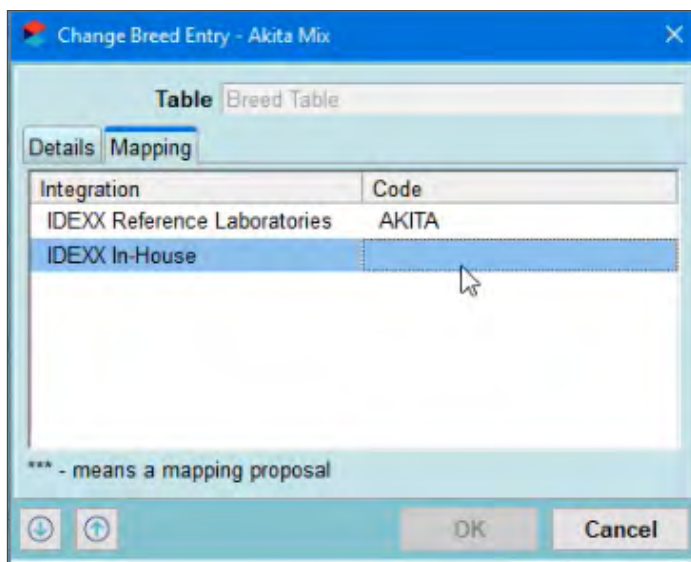


## Select the Device

When there are multiple IVLS machines, a prompt will appear for the user to select the desired device for sending the lab request. Avimark will remember the last IVLS selected on the computer but a different device can be selected from the drop-down list.

**Note:** This prompt will **not** appear if there is only one IVLS.


1. In the Patient area, **right-click > Create Requisition.**
2. Select the IVLS device from the drop-down list.

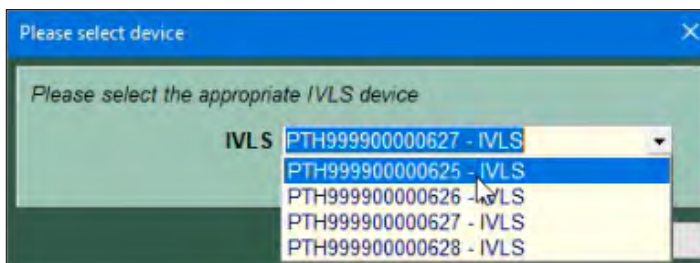

3. Proceed with placing the order.

## Place the Order


After the IVLS device is selected, the screen will appear for the user to click the button to place the order.

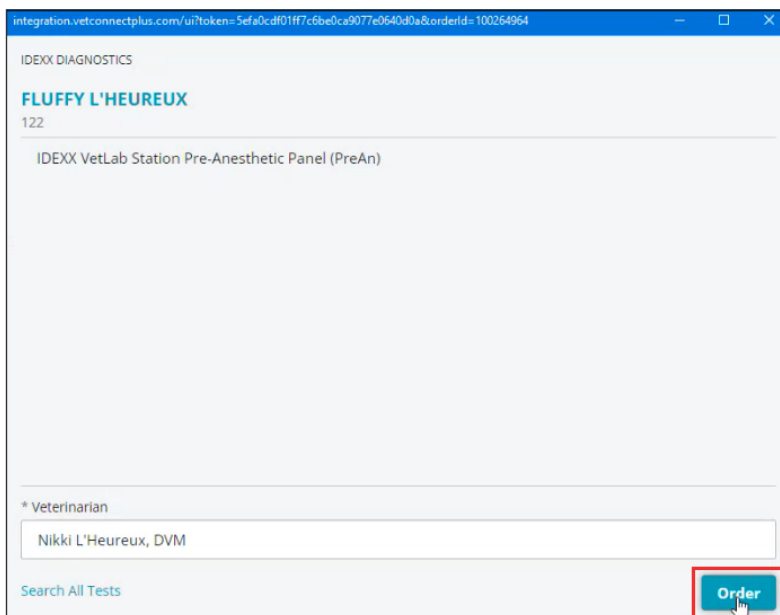

## Notifications

IDEXX In-House lab results, when complete, will appear under Notifications.

1. Right-click on a completed result.
2. Click on **Import Results**.

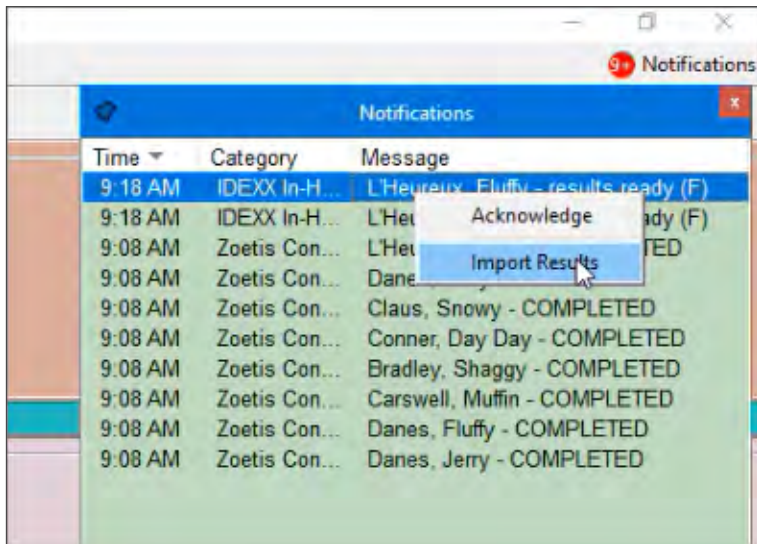

3. Requisition Management will load the IDEXX In-House labs and download the results.

